

# *Pediatric Crib* *For General and Intensive Care*

MODEL ICC 411



**SAVION**  
ADVANCED PATIENT CARE®

## Features - Benefits

- unique self-closing and automatic locking oval shape sidebars
- siderails operated by locking pneumatic gas springs with any number of intermediate stop positions
- hands-free foot-pedal operated siderail adjustment – cannot be accidentally released by the child
- two-way tilting steel mesh platform, infinitely adjustable by pneumatic gas springs – for trendelenburg, reverse trendelenburg and adjustment of platform height
- full-length underbed ABS cover with built-in storage tray
- smooth-rolling two-tone twin-wheel castors – all individually braked
- Mattress – 4" (10cm) pressure sore reducing poly-foam, with waterproof, anti-bacterial, anti-static, fire-retardant stretch cover
- Pkd delivery available for on-site assembly
- optional: electrically operated platform tilt and height adjustment, with push-button handset control

## General

The ICC411 pediatric crib/cot sets new standards for pediatric care. All-round access is available when needed, whilst sides can be raised


effortlessly through the use of the pedal operated gas-spring assisted mechanism and therefore require no lifting effort by attendants and nurses. At the same time, any number of intermediate stops are available for every type of situation and requirement. At every intermediate height, the locking is total, cannot be opened by the child inside the bed and is perfectly safe.

Sides cannot drop suddenly even if the nurse lets go of the rail before it is locked.

If trendelenburg or reverse trendelenburg is needed – this can be performed at either end of the crib and allows the platform to be set at an any number of height positions to suit every



any number of intermediate siderail stop positions

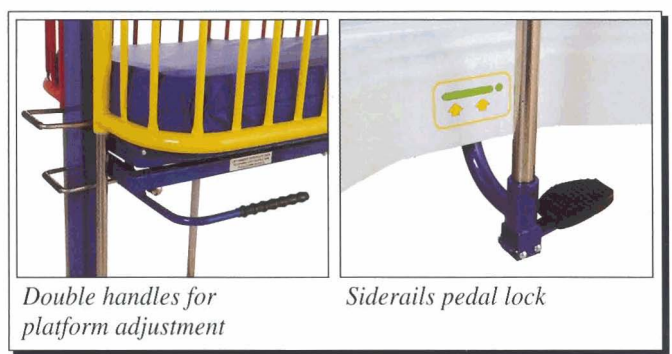

Double handles for platform adjustment

Siderails pedal lock

medical situation, child age or nursing requirement. Platform is manufactured from epoxy coated steel mesh ensuring maximum below-mattress ventilation. Below the platform is a stylish full-length ABS cowling, protecting the underbed mechanism and providing a large dished storage tray for diapers, other patient effects or equipment.

## Construction

The entire chassis is manufactured from welded tubular steel and can be delivered in knock-down form for reassembly without the need for any special tools or equipment.

Castors are 5" twin-wheel rust-free plastic, which combine strength and enhanced mobility. All 4 castors have an easy to operate foot-operated double locking mechanism.

## Finish

Cribs are available in a wide range of attractive colour finishes. Highest grade epoxy paint, electrostatically sprayed and baked at high temperature, provides a rich and pleasant finish, resistant to chipping and scratching.

## Accessories

An IV drip pole is fitted as standard in one of the 2 accessory wells.

Pump poles as well as a range of other accessories are also available.

## Dimensions

Mattress Platform			Overall		Trendelenburg and reverse trendelenburg	
92/138cm 36/54"	54/73cm 21.5/29"	100cm - 81cm 39.5" - 32"	108/154cm 42.5/60.5"	64/83cm 25/32.5"	10°	10°

Encore Medical Technologies, Inc.  
12804 S. Spring Street, Los Angeles, CA 90061  
P: 800-558-8636, 310-538-9138, F: 310-538-9177  
[www.encoremedical.com](http://www.encoremedical.com)  
Email: [info@encoremedical.com](mailto:info@encoremedical.com)

Specifications are subject to change from time to time

